


Electrosurgical Grounding Pads

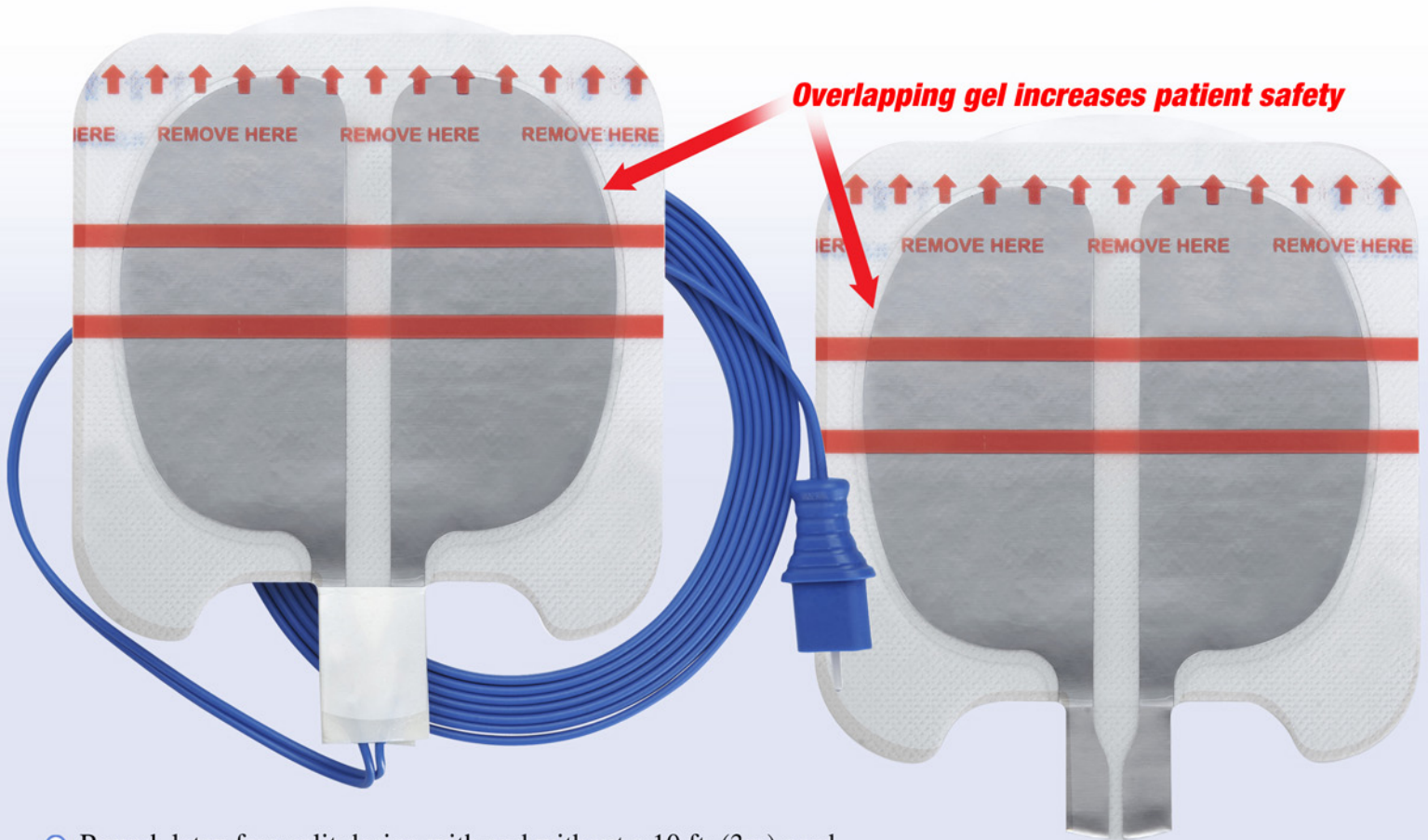

- Round, latex free split design with and without a 10 ft. (3m) cord
- Overlapping gel-line border reduces risk of incident as aluminum edge exposure is eliminated
- Coated non-woven material allows for greater heat dissipation
- Designed for compatibility with all industry standard generators that meet the ANSI/AAMI standard HF 18:2001 for Electrosurgical Devices
- Compatible with generators designed with a contact quality monitor, deactivating if the patient impedance is reduced at the grounding pad
- Compact size: 14 in² (90 cm²) – Enables versatile placement in any direction
- For use on patients over 55 pounds (25 kg)

Electrosurgical Grounding Pads

Product #	Product Description
20240	Split Grounding Pad with 10 ft. Cord
20230	Split Grounding Pad without Cord



Pittsburgh, PA U.S.A • Tel: 724-337-5500 Tel: 800-963-8776
Fax: 724-337-0555 Fax: 800-963-6553

www.xodusmedical.com • email: info@xodusmedical.com