


Electrosurgical Grounding Pads



- Round, latex free split design with and without a 10 ft. (3m) cord
- Overlapping gel-line border reduces risk of incident as aluminum edge exposure is limited
- Coated non-woven material allows for greater heat dissipation
- Compatible with generators designed with a contact quality monitor, deactivating if patient impedance is reduced at the grounding pad
- For use on patients over 55 pounds (25 kg)

Product #	Product	Case Qty
20240	Split Grounding Pad with Cord	50
20230	Split Grounding Pad without Cord	50

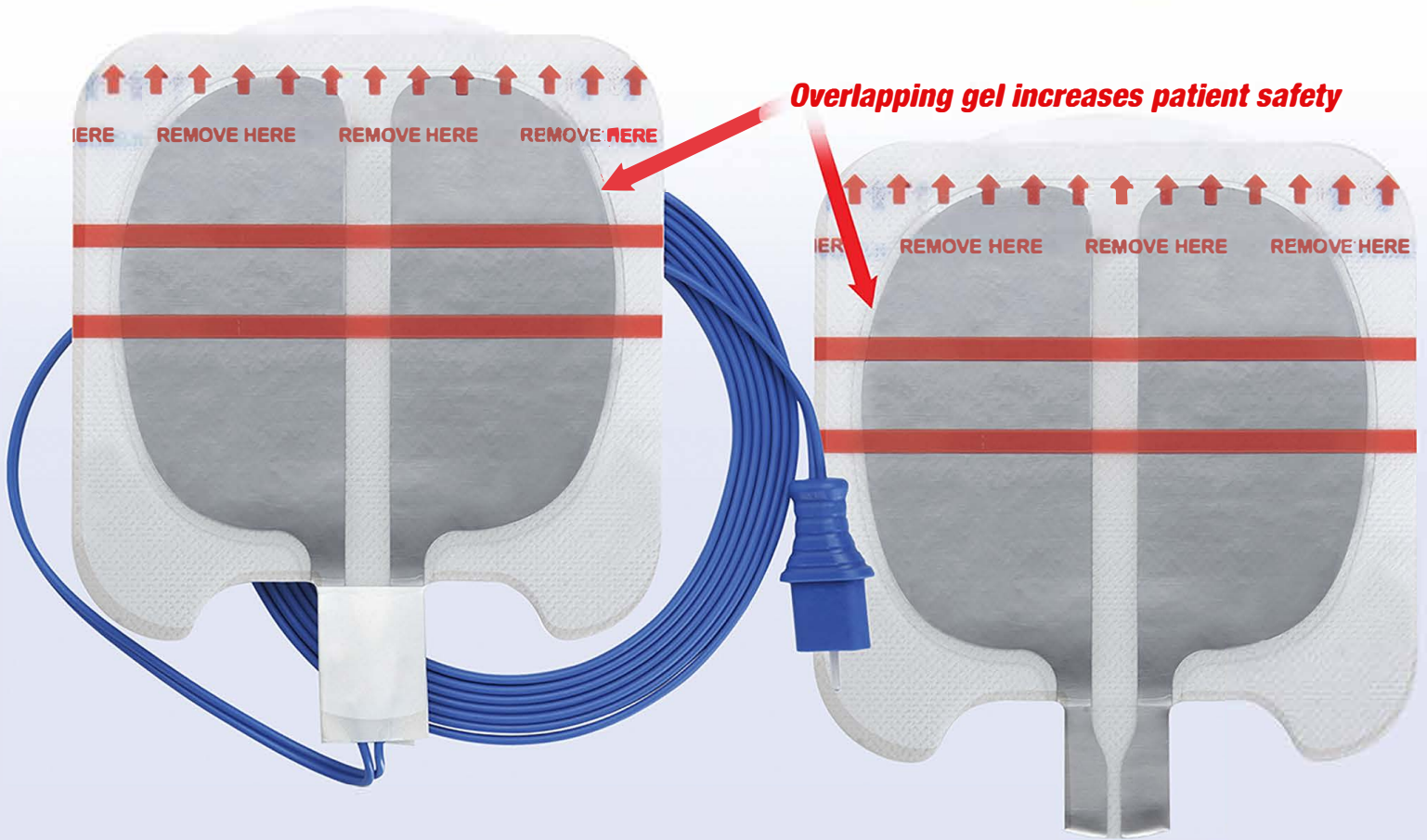

ISO 13485 CERTIFIED
FDA REGISTERED

Pittsburgh, PA USA • Tel: 724-337-5500
Fax: 724-337-0555

xodusmedical.com
info@xodusmedical.com



LITBRES003 R3 11/26/18